

CPAP VT 80D *Aqua*



iPR (Intelligent Pressure Relief)

This feature allows users relief when exhaling against the pressure. There are three iPR levels to meet your individual needs.

Auto Start/ Stop

This feature will turn the device **On** for you have the mask on without the need to press any buttons and will also turn the device **Off** once you take your mask Off.


Technical Specification

Modes	CPAP, Auto CPAP	
Operating/ Monitoring Range	4-20 CmH2O	
Sound Pressure Level	<28 dBA with uncertainty of 2 dBA	
Leak Compensation	Auto	
Ramp Time	0-60 mins	
iPR	On/ Off, 1-3 levels	
Humidification Level	Off, 1-6 level	
Display	3.5 in LCD	
Alarms	Leak Alert	
Data Management	1/7/31 days, 6 months, 1 year	
Dimensions (LxWxH)	282x172.5x141.5 mm (Device with Humidifier)	
Weight	1900 g (Device with Humidifier)	
Power Adapter	Input Range : 100-240V, 50/60 H ² , 1.5A Max, Over voltage protection Auto Controlled Temp variation	
Battery Backup	External 2-8 Hrs (Optional)	
Operating Temperature/ Humidity/ Altitude	+5 °c to +35 °c/ 10-95% relative humidity, non-condensing/sea level to 2,591 m, air pressure range 70 kPa to 105 kPa	
Storage and Transport Temperature/Humidity	-20 °c to +60 °c / 5-95% relative humidity, non-condensing	
Supplement Oxygen	Maximum Flow: 4L/min (CPAP)	
Standard Air Filter	Reusable Cotton Filter	
Tubing	Heating Tube (Optional)	
Air Outlet	22 mm Air Outlet (STD)	
Warranty	2 Years	
Conforms to Standards	IEC 60601-1-1 IEC 60601-1-2 IEC 80601-2-12 CDSCO	Class 2 (double insulated, Type BF, Ingress Protection IP22)

*Term & Conditions Apply

Hitech Turbine Excellent Performance

- Idler turbine & Quick Response
- Low-noise operation and deep tranquility

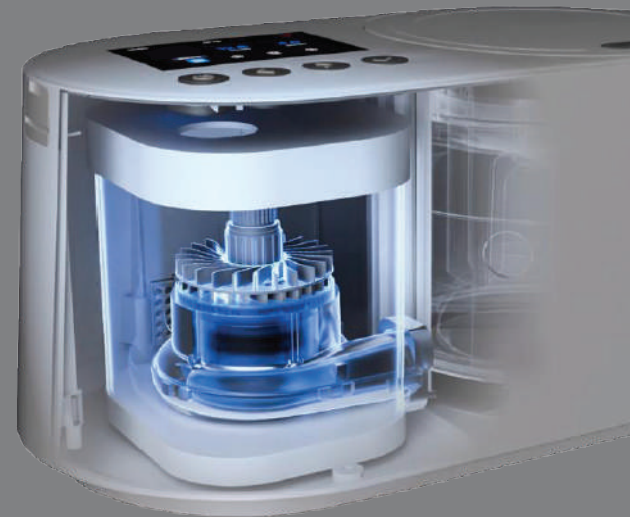

Integrated Automatic Humidification

Deck Mount unique Intelligent humidifier can apperceive the environment temperature & humidity and maintain the humidification within certain level so that there will not be condensation in the tube and mask.

Product Code

Respirare Auto CPAP with Humidifier

VT 80 D

Note: Technical specification may change without notice. Always refer to your device user manual.

Manufactured by:

Deck Mount Electronics Pvt. Ltd 260, Phase IV, Udyog Vihar, Sector 18, Gurugram, Haryana 122015

1800 30 92 499

Made in India

wecare@deckmount.in

www.deckmount.in

DECK MOUNT
Saving millions of lives

 **respirare**

Silent Sleep. Smarter Therapy



RESPIRARE VT 80D *Aqua*

Experience a new era of sleep therapy with Respirare. Our advanced CPAP device delivers exceptional quietness, combined with cutting-edge monitoring and comprehensive reporting.