

Intelligent Breathing Auto CPAP Therapy

Breathing Innovation, Personalized for You!



RESPIRARE VT 10D


A top-quality CPAP machine delivers reliable, adjustable pressure levels with advanced comfort features and quiet operation for effective respiratory support.

CPAP VT 10D



Alarm Function

The alarm level meets national standards, sound alarm + graphic alarm, timely remind the cause of the alarm, for the safety of users.

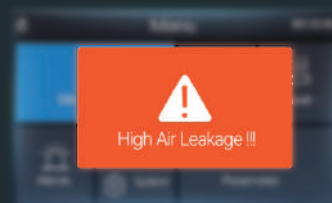

Technical Specification

Modes	CPAP, APAP, iAPAP	
Operating/ Monitoring Range	4-20 CmH2O	
Sound Pressure Level	<28 dBA with uncertainty of 2 dBA	
Leak Compensation	Auto	
Auto On,Off/ SD card	On/ Off	
Ramp Time/ i Ramp	0-60 mins/ On/ Off	
B-Flex	0-3 levels	
Preheat/ Preheat time	On/ Off / 10-30 mins	
Humidification Level	0-5 level	
Display	2.8 in LCD	
Boost Sensitivity	Standard Sensitivity Soft	
Alarm	High Leakage Apnea	
Dimensions (LxWxH)/ Weight	250x181x114 mm/ 1340 g	
Data Management	1/7/31 days, 6 months, 1 year of detailed data.	
Power Adapter	Input Range : 100-240V, 50/60 Hz, 1.8 A max (typical power consumption)	
Operating Temperature/ Humidity/ Altitude	+5 °c to +35 °c/ 10-95% relative humidity, non-condensing/sea level to 2,591 m, air pressure range 700 to 1060 hpA	
Storage and Transport Temperature/Humidity	-20 °c to +60 °c / 5-95% relative humidity, non-condensing	
Supplement Oxygen	Maximum Flow: 4L/min (CPAP)	
Standard Air Filter	Cotton Filter	
Tubing/ Air Outlet	Standard (19 mm)/ 22 mm Air Outlet (STD)/ Heating tube (Optional)	
Sound Pressure Level	<30 dBA with uncertainty of 2 dBA	
Warranty	2 Years (Repair)	
Conforms to Standards	IEC 60601-1-1 IEC 60601-1-2 IEC 80601-2-12 CDSCO	Class 2 (double insulated, Type BF, Ingress Protection IP22)

*Term & Conditions Apply

Key Features

B-Flex

The B-Flex is a comfort feature that will lower the pressure when you exhale and slowly increases the pressure when you inhale and will be available in Auto CPAP mode.

Auto On/Off

This feature will turn the device **On** for you have the mask on without the need to press any buttons and will also turn the device **Off** once you take your mask Off.

Integrated Humidification

Heated Humidification makes therapy more comfortable and can reduce side effects like a dry nose and throat.

Product Code

Respirare Auto CPAP with Humidifier VT 10 D

Note: Technical specification may change without notice. Always refer to your device user manual.

Deck Mount Electronics Pvt. Ltd

260, Phase IV, Udyog Vihar, Sector 18, Gurugram, Haryana 122015

1800 30 92 499

Made in India

wecare@deckmount.in

www.deckmount.in